




## FLEXIBLE STRENGTH. PROVEN SUCCESS.

Demonstrated long-term durable clinical outcomes in complex aortoiliac occlusive disease (AIOD) treatment through **5 years**.

**Sustained clinical effectiveness through 5 years:**

- 89.5% primary patency and 96.1% primary assisted patency per lesion<sup>1</sup>
- 89.1% freedom from target lesion revascularization (fTLR) per subject<sup>1</sup>

This physician-initiated study enrolled 59 patients from 3 participating centers with patients followed out to **5 years**.

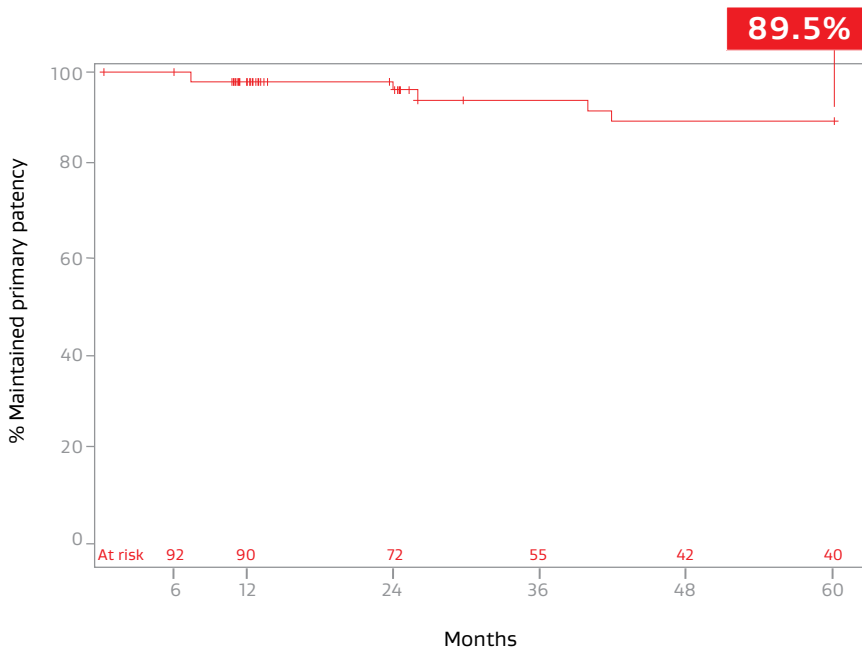

Kaplan-Meier graph of primary patency with number of lesions at risk.

# DURABLE PATIENT BENEFIT VERSUS BASELINE AT 5 YEARS<sup>1</sup>

**100%** of evaluated (n=28) patients  
improved  $\geq 1$  Rutherford category  
from baseline<sup>\*,1</sup>

**.15** improvement in mean resting ankle-brachial  
index (ABI) ( $P < .001$ , .95 mean ABI)<sup>†,‡</sup>

3x improvement in median WIQ score sustained  
beyond 5 years in long-term follow-up cohort<sup>1</sup>

	Pre-procedure	3 years	5 years <sup>‡</sup>
Walking distance	7 (N = 59)	25 (N = 39)	21 (N = 27)
Walking speed	3 (N = 59)	10 (N = 39)	9 (N = 27)
Stair climbing	3 (N = 50)	11 (N = 38)	9 (N = 27)

\* 59 subjects participated and 28 were available through the end of the study at 5-year follow-up.

† ( $P < .001$ ) Statistically significant change from pre-procedure.

‡ Median follow-up of 6.6 years.



## Reference

1. Holden A, Takele E, Hill A, et al. Long-Term Follow-up of Subjects With Iliac Occlusive Disease Treated With the Viabahn VBX Balloon-Expandable Endoprosthesis. *Journal of Endovascular Therapy* 2023;0(0). Doi:10.1177/15266028231165723

 Consult Instructions  
for Use  
[eifu.goremedical.com](http://eifu.goremedical.com)

Refer to *Instructions for Use* at [eifu.goremedical.com](http://eifu.goremedical.com) for a complete description of all applicable indications, warnings, precautions and contraindications for the markets where this product is available. <sup>Rx only</sup>

Products listed may not be available in all markets.

GORE, *Together, improving life*, VBX, VBX FLEX, VIABAHN and designs are trademarks of W. L. Gore & Associates.  
© 2024 W. L. Gore & Associates GmbH 24PL2041-EN01 MAY 2024

W. L. Gore & Associates, Inc.  
[goremedical.com](http://goremedical.com)

Asia Pacific +65 67332882 Australia / New Zealand 1800 680 424 Europe 00800 6334 4673  
United States Flagstaff, AZ 86004 800 437 8181 928 779 2771

