




FLEXIBLE STRENGTH. PROVEN SUCCESS.

Demonstrated long-term durable clinical outcomes in complex aortoiliac occlusive disease (AIOD) treatment through **five years**.

Sustained clinical effectiveness through five years:

- 89.5% primary patency and 96.1% primary assisted patency per lesion¹
- 89.1% freedom from target lesion revascularization (fTLR) per subject¹

This physician-initiated study enrolled 59 patients from three participating centers with patients followed out to **five years** and beyond.

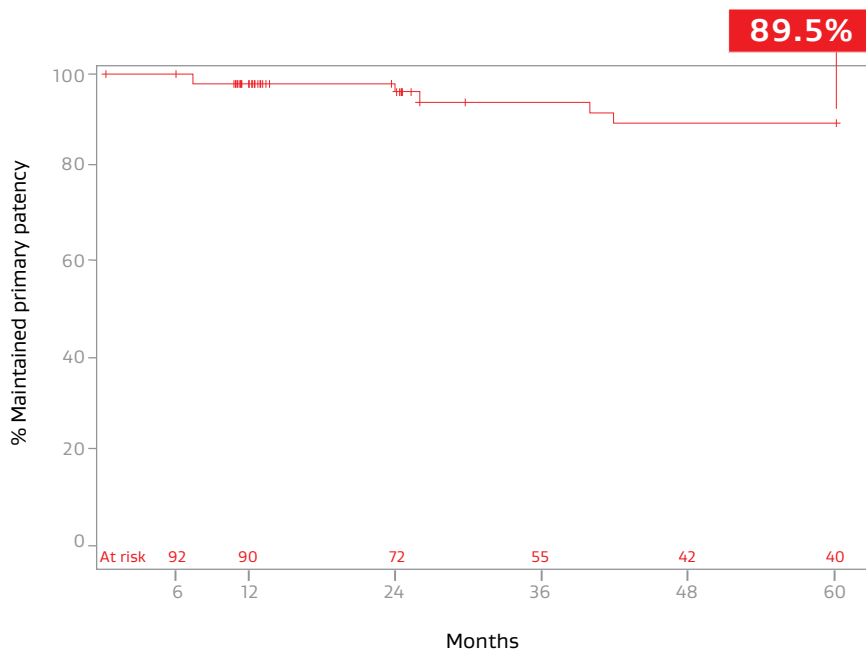

Kaplan Meier graph of primary patency with number of lesions at risk

DURABLE PATIENT BENEFIT VERSUS BASELINE BEYOND FIVE YEARS¹

89% of patients improved ≥ 1
Rutherford category from baseline¹

.15 improvement in mean resting ankle-brachial
index (ABI) ($P < .001$, .95 mean ABI)^{*,†}

3x improvement in median WIQ measures sustained
beyond five years in long-term follow-up cohort¹

	Pre-procedure	3 years	5 years [†]
Walking distance	7 (N = 59)	25 (N = 39)	21 (N = 27)
Walking speed	3 (N = 59)	10 (N = 39)	9 (N = 27)
Stair climbing	3 (N = 50)	11 (N = 38)	9 (N = 27)

* ($P < .001$) Statistically significant change from pre-procedure.

† Median follow-up of 6.6 years.



Reference

1. Holden A, Takele E, Hill A, et al. Long-term follow-up of subjects with aortoiliac occlusive disease treated with the VIABAHN VBX balloon expandable endoprosthesis. Presented at Cardiovascular & Interventional Radiological Society of Europe (CIRSE) 2022; September 10-14, 2022; Barcelona, Spain. *CardioVascular & Interventional Radiology* 2022;45:Supplement:S147. Abstract 301.1.

 Consult Instructions
for Use
eifu.goremedical.com

Refer to *Instructions for Use* at eifu.goremedical.com for a complete description of all applicable indications, warnings, precautions and contraindications for the markets where this product is available. ^{Rx only}

Products listed may not be available in all markets.

GORE, *Together, improving life*, VBX, VBX FLEX, VIABAHN and designs are trademarks of W. L. Gore & Associates.
© 2022 W. L. Gore & Associates GmbH 22757249-EN NOVEMBER 2022

W. L. Gore & Associates, Inc.
goremedical.com

Asia Pacific +65 67332882 Australia / New Zealand 1800 680 424 Europe 00800 6334 4673
United States Flagstaff, AZ 86003 800 437 8181 928 779 2771

