




Together, improving life

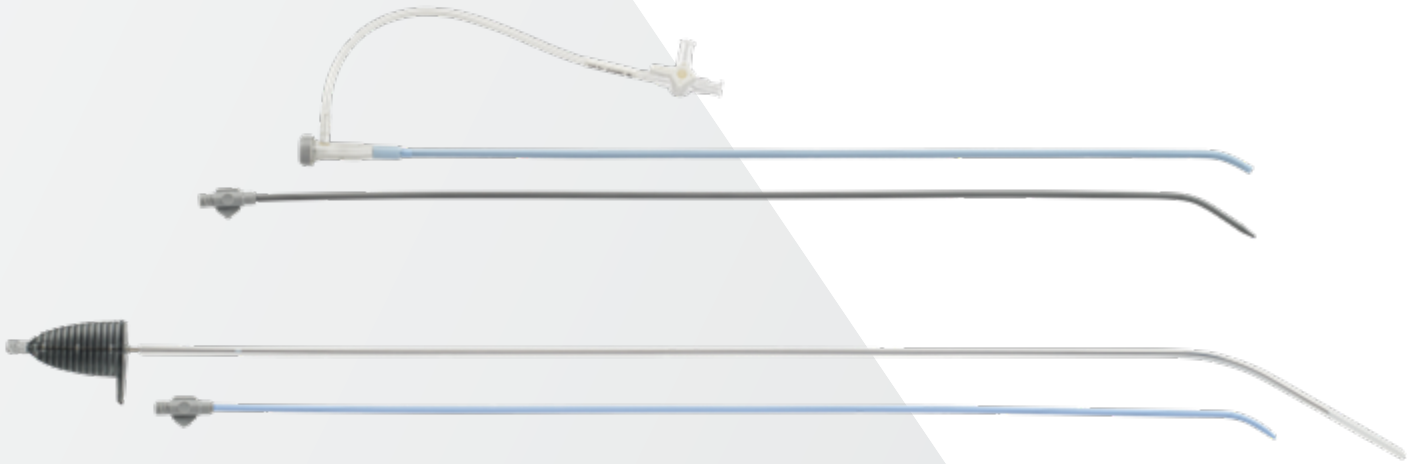


GORE TIPS Set

DESIGNED TO DELIVER




Advanced accessories
exclusively tailored for
TIPS procedures.

GORE TIPS Needle

Ergonomic handle enhances grip and needle control.

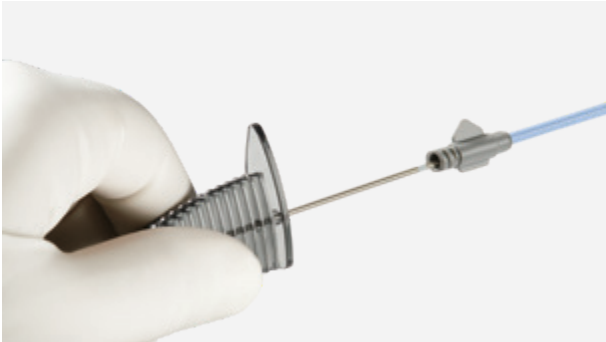

Direction indicators verify component alignment during device preparation and initial insertion.

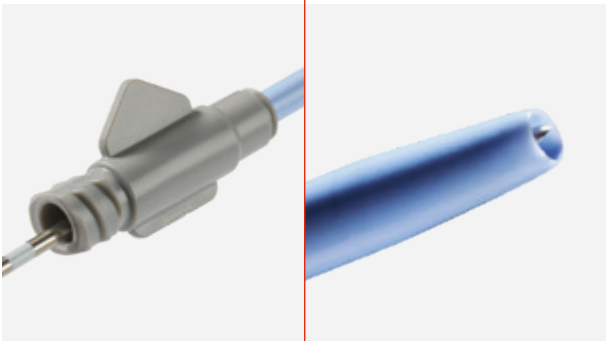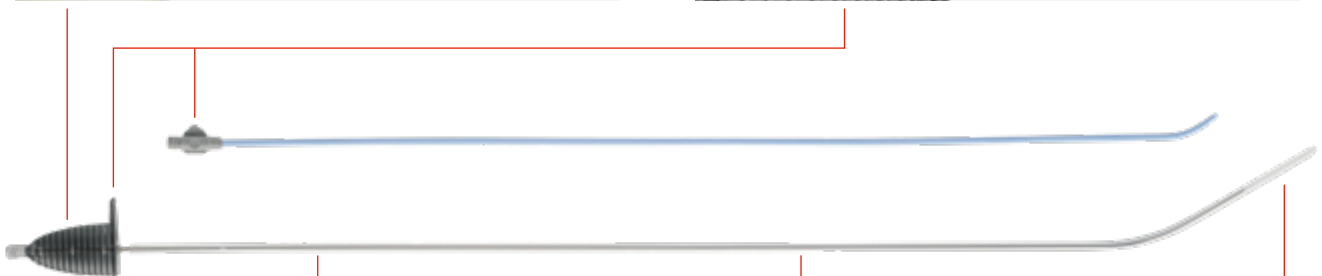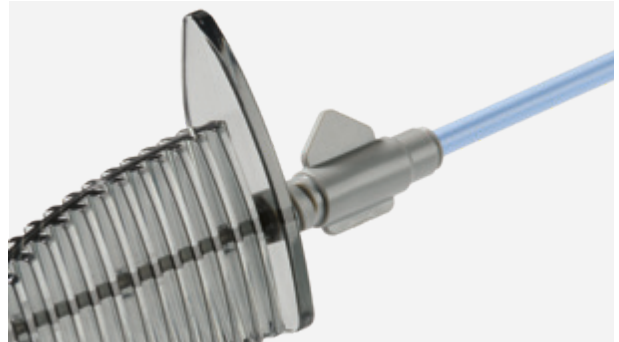

Status band confirms needle tip protection.



Stiff 16 Gauge needle provides responsive needle torque control.

The GORE TIPS Needle features a status band that when visible confirms that the needle tip is covered.




Sharp needle tip aids precise needle puncture.

GORE® TIPS Sheath


Interlocking feature secures dilator hub and 10 Fr introducer sheath during advancement.



Transparent tapered hub provides clear view of access sleeve and device delivery catheter.

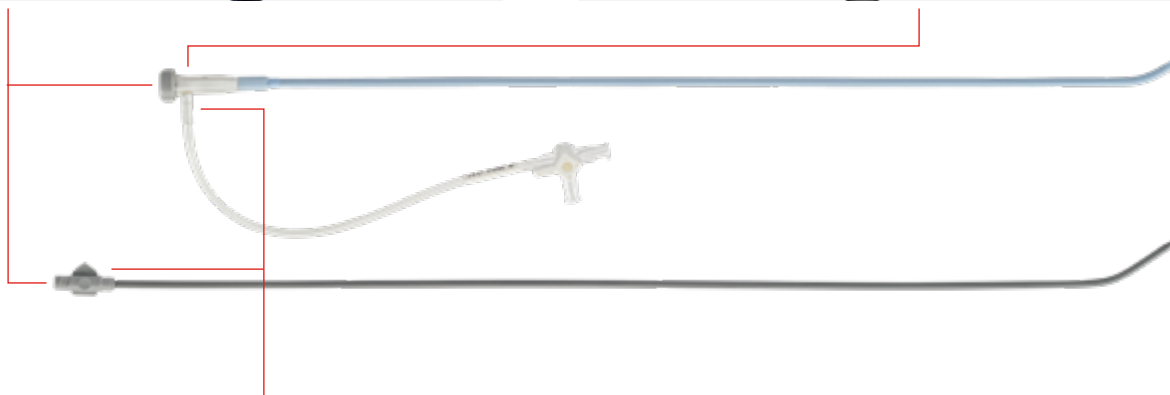

Direction indicators verify component alignment during device preparation and initial insertion.

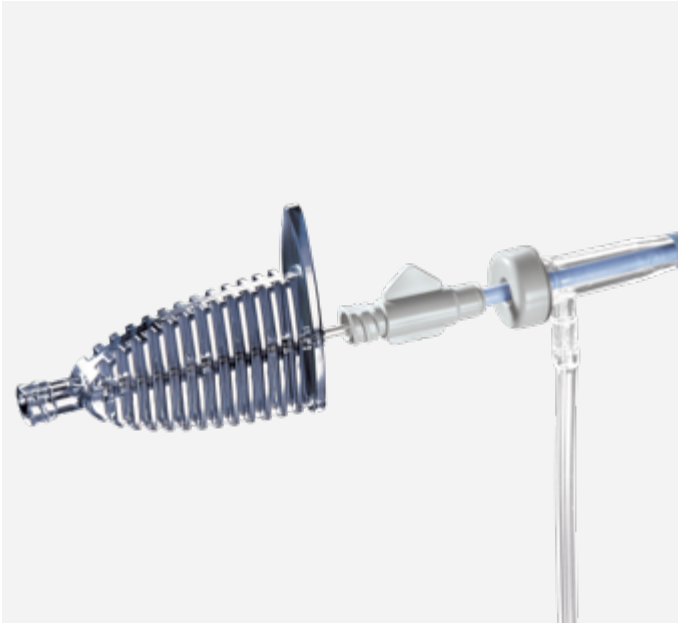

The GORE TIPS Set consists of the GORE® TIPS Sheath and the GORE TIPS Needle.



The GORE TIPS Needle consists of a high-strength, stiff 16 Gauge stainless steel hollow shaft needle and 10 Fr guiding catheter.



The GORE® TIPS Sheath consists of a flexible 10 Fr introducer sheath with a hemostatic valve and a 10 Fr dilator.

Configurations

GORE TIPS Set

Catalogue number	GORE TIPS Set component	Effective length (cm)	Accepts guidewire diameter (in)
TSET1016	16 Gauge needle	56	≤ .035
	10 Fr Introducer sheath	40	≤ .035
	10 Fr Dilator	47	≤ .035
	10 Fr Guiding catheter	49	≤ .035

 Consult Instructions for Use
eifu.goremedical.com

Refer to *Instructions for Use* at eifu.goremedical.com for a complete description of all applicable indications, warnings, precautions and contraindications for the markets where this product is available. Rx Only

Products listed may not be available in all markets.

GORE, *Together, improving life*, and designs are trademarks of W. L. Gore & Associates.
 © 2022 W. L. Gore & Associates GmbH 22687257-EN AUGUST 2022

W. L. Gore & Associates, Inc.
goremedical.com

Asia Pacific +65 6733 2882 Australia/New Zealand 1800 680 424 Europe 00800 6334 4673
 United States Flagstaff, AZ 86004 800 437 8181 928 779 2771

